

ZytoDot® SPEC MDM2 Probe




Background

The ZytoDot® SPEC MDM2 Probe (PD9) is intended to be used for the qualitative detection of amplifications involving the human MDM2 gene in formalin-fixed, paraffin-embedded specimens, such as atypical lipomatous tumor/well-differentiated liposarcoma (ALT/WDLPS) and dedifferentiated liposarcoma (DDLPS), by chromogenic *in situ* hybridization (CISH). The probe is intended to be used in combination with the ZytoDot® CISH Implementation Kit (Prod. No. C-3018-40). The product is intended for professional use only. All tests using the product should be performed in a certified, licensed anatomic pathology laboratory under the supervision of a pathologist/human geneticist by qualified personnel. The probe is intended to be used as an aid to the differential diagnosis of ALT/WDLPS and DDLPS and therapeutic measures should not be initiated based on the test result alone.


Probe Description

The ZytoDot® SPEC MDM2 Probe is composed of:

- Digoxigenin-labeled polynucleotides (~1.8 ng/μl), which target sequences mapping in 12q15** (chr12:69,190,708-69,430,185) harboring the MDM2 gene region.
- Formamide based hybridization buffer

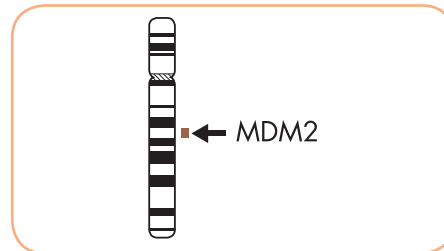


Ideogram of chromosome 12 indicating the hybridization locations.

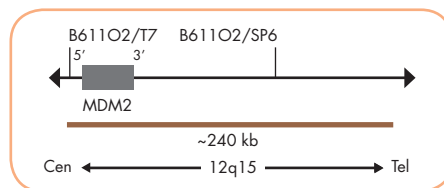

SPEC MDM2 Probe map (not to scale).

Results

In normal cells, two distinct dot-shaped signals per nucleus will be observed. Nuclei with amplification of the MDM2 gene locus or polysomy of chromosome 12 will show multiple dots or large signal clusters.

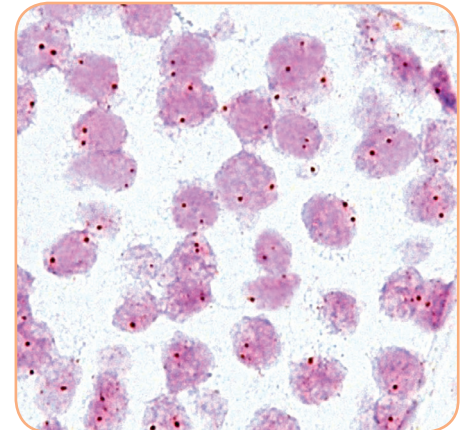

Normal nuclei each with two MDM2 signals.

Prod. No.	Product	Label	Tests* (Volume)
C-3012-400	ZytoDot SPEC MDM2 Probe	DIG	40 (400 μl)
Related Products			
C-3018-40	ZytoDot CISH Implementation Kit		40
<small>Incl. Heat Pretreatment Solution EDTA, 500 ml; Pepsin Solution, 4 ml; Wash Buffer SSC, 560 ml; PBS/Tween, good for 2000 ml; Blocking Solution, 4 ml; Mouse-anti-DIG, 4 ml; Anti-Mouse-HRP-Polymer, 4 ml; DAB Solution A, 0.3 ml; DAB Solution B, 10 ml; Mayer's Hematoxylin Solution, 20 ml; Mounting Solution (alcoholic), 4 ml</small>			

* Using 10 μl probe solution per test. labeled products are only available in certain countries. All other countries research use only! Please contact your local dealer for more information.

**According to Human Genome Assembly GRCh37/hg19